

UniGo Microfluidic Pump

Accurate and Pulse-Free Flow for Microfluidic Applications

🎯 SUPERIOR PERFORMANCE

UniGo microfluidic pump is a precision pressure pump for a variety of microfluidic applications, where accurate and stable flow rate delivery is required.

🏗️ MODULAR

Up to 4 UniGo microfluidic pumps racked together allowing simultaneous control and independent programming of each pump's flow profile.

😊 EASY TO USE

Control all the features of the pump with the SmartFlo app running on iPad mini or PC.



☰ MAIN BENEFITS



FAST

Fast response time and settling time



PROGRAMMABLE

Program your own flow profile with any combination of constant, ramp, sine and step.



PULSE FREE

Achieve precise and pulse-less flow.



EASY INTEGRATION

Integrate the pump control on your own software using LabVIEW, C++, Python, etc.



VERSATILE

Combine it with ExiGo pumps to get the best features of both pressure and syringe pumps in one microfluidic set-up.




CLOSED-LOOP FEEDBACK

Integrates a plug-and-play flow sensor for active feedback and precise flow control.



How does it work?

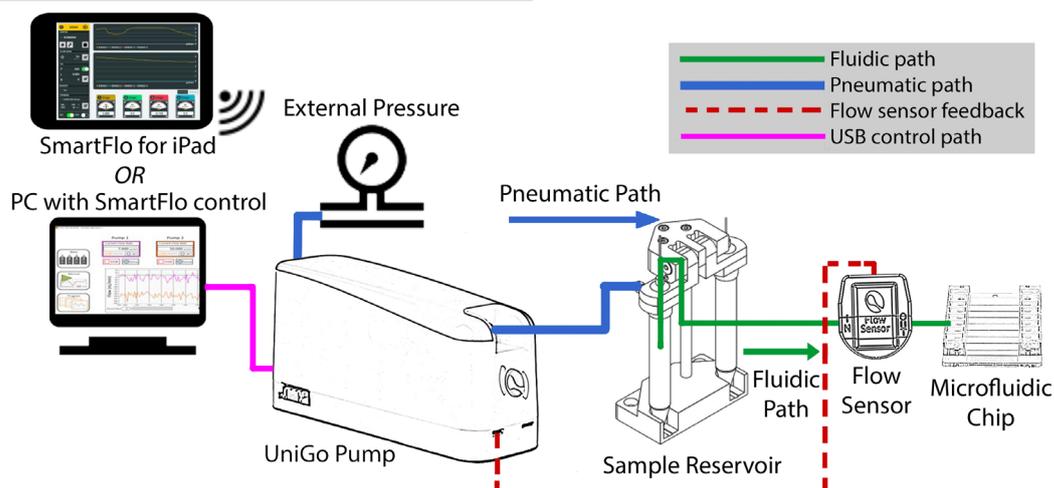

TECHNICAL SPECIFICATIONS

Input Pressure	Up to 10 bar or 145PSI or vacuum -0.9bar ⁽¹⁾
Volumetric Flow Rates	1µL/min to 1mL/min
Flow Direction	Unidirectional (Push or Pull)
Pulse-free flow rate for microfluidic applications	1µL/min - 80µL/min (\pm 120nL/min) ⁽²⁾ 10µL/min - 1mL/min (\pm 2µL/min) ⁽³⁾
Flow feedback options	Pressure-flow combined or Pressure only
Fast response	<50ms
Software Control	SmartFlo App running on iPad mini or PC
Dimensions	225mm (L) x 69mm (W) x 122mm (H)
Weight	~1.3kg
Power Requirements	110/220V - 50/60Hz 60W

(1) External compressor or vacuum required.

(2) Using a 80µL/min sensor (FS80.0)

(3) Using a 1mL/min sensor (FS-1000).

APPLICATIONS

- ✓ Droplet Generation
- ✓ Organ-on-chip studies
- ✓ VOC experimental set-ups
- ✓ Shear based cell rolling, adhesion and transmigration assays on ligand-coated surfaces or endothelial cells
- ✓ Cell and particle manipulation
- ✓ Precise multichannel mixing
- ✓ Bacterial biofilm formation on solid surfaces

What's included?

	INCLUDED	OPTIONAL
UniGo pump	✓	
Flow Sensor ⁽⁴⁾	✓	
iPad Mini ⁽⁵⁾		✓
Sample Reservoir		✓
Fluidic Manifold		✓

(4) More information available on www.wearecellix.com/flowsensors

(5) SmartFlo is controlled via iPad or PC. Please contact Cellix for more information.